



## Research article

# Improving *gonorrhoea* molecular diagnostics: Genome mining-based identification of identical multi-repeat sequences (IMRS) in *Neisseria gonorrhoeae*

Clement Shiluli<sup>a</sup>, Shwetha Kamath<sup>b</sup>, Bernard N. Kanoi<sup>a</sup>, Racheal Kimani<sup>a</sup>, Michael Maina<sup>a</sup>, Harrison Waweru<sup>a</sup>, Moses Kamita<sup>a</sup>, Ibrahim Ndirangu<sup>a</sup>, Hussein M. Abkallo<sup>d</sup>, Bernard Oduor<sup>d</sup>, Nicole Pamme<sup>e</sup>, Joshua Dupaty<sup>c</sup>, Catherine M. Klapperich<sup>c</sup>, Srinivasa Raju Lolabattu<sup>b,\*\*</sup>, Jesse Gitaka<sup>a,\*</sup>

<sup>a</sup> Centre for Research in Infectious Diseases, College of Graduate Studies and Research, Mount Kenya University, Thika, Kenya

<sup>b</sup> Division of Research and Development, Jigsaw Bio Solutions Private Limited, Bangalore, India

<sup>c</sup> Department of Biomedical Engineering, Boston University, Boston, MA, USA

<sup>d</sup> Animal and Human Health Program, International Livestock Research Institute, Nairobi, Kenya

<sup>e</sup> Department of Materials and Environmental Chemistry, Stockholm University, Sweden



## ARTICLE INFO

## Keywords:

Identical multirepeat sequences  
*Neisseria gonorrhoeae*  
Isothermal

## ABSTRACT

**Background:** Curable sexually transmitted infections (STIs), such as *Neisseria gonorrhoeae* (*N. gonorrhoeae*), are a major cause of poor pregnancy outcomes. The infection is often asymptomatic in pregnant women, and a syndrome-based approach of testing leads to a missed diagnosis. Culture followed by microscopy is inadequate and time-consuming. The gold standard nucleic acid amplification tests require advanced infrastructure settings, whereas point-of-care tests are limited to immunoassays with sensitivities and specificities insufficient to accurately diagnose asymptomatic cases. This necessitates the development and validation of assays that are fit for purpose.

**Methods:** We identified new diagnostic target biomarker regions for *N. gonorrhoeae* using an algorithm for genome mining of identical multi-repeat sequences (IMRS). These were then developed as DNA amplification primers to design better diagnostic assays. To test the primer pair, genomic DNA was 10-fold serially diluted (100 pg/μL to  $1 \times 10^{-3}$  pg/μL) and used as DNA template for PCR reactions. The gold standard PCR using 16S rRNA primers was also run as a comparative test, and both assay products were resolved on 1% agarose gel.

**Results:** Our newly developed *N. gonorrhoeae* IMRS-PCR assay had an analytical sensitivity of 6 fg/μL representing better sensitivity than the 16S rRNA PCR assay with an analytical sensitivity of 4.3096 pg/μL. The assay was also successfully validated using clinical urethral swab samples. We further advanced this technique by developing an isothermal IMRS, which was both reliable and sensitive for detecting cultured *N. gonorrhoeae* isolates at a concentration of 38 ng/μL. Combining isothermal IMRS with a low-cost lateral flow assay, we were able to detect *N. gonorrhoeae* amplicons at a starting concentration of 100 pg/μL.

\* Corresponding author.

\*\* Corresponding author.

E-mail addresses: [raju@jigsawbio.com](mailto:raju@jigsawbio.com) (S.R. Lolabattu), [jgitaka@mku.ac.ke](mailto:jgitaka@mku.ac.ke) (J. Gitaka).

<https://doi.org/10.1016/j.heliyon.2024.e27344>

Received 24 August 2023; Received in revised form 19 February 2024; Accepted 28 February 2024

Available online 5 March 2024

2405-8440/© 2024 Published by Elsevier Ltd.

This is an open access article under the CC BY-NC-ND license

(<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

**Conclusion:** Therefore, there is a potential to implement this concept within miniaturized, isothermal, microfluidic platforms, and laboratory-on-a-chip diagnostic devices for highly reliable point-of-care testing.

## 1. Introduction

Gonorrhoea caused by *Neisseria gonorrhoeae* (*N. gonorrhoeae*) accounted for an estimated 82.4 million cases globally in 2020 [1]. Problems associated with *N. gonorrhoeae* infections include ectopic pregnancy, pelvic inflammatory disease, and infertility [2].

Microscopy provides an affordable diagnostic option for symptomatic infections, but requires trained personnel and has low sensitivity in asymptomatic individuals [3,4]. In men, *N. gonorrhoeae* microscopy is 84%–95% sensitive compared with culture, and it is more sensitive in symptomatic (90%–95%) than in asymptomatic (50%–75%) patients [5]. In female patients, microscopic observation of cervical smears identifies approximately 50% of *N. gonorrhoeae* infections and is therefore not recommended for routine testing [6–8].

Antigen-based tests detect *N. gonorrhoeae* using antibodies immobilized on nitrocellulose strips [9]. However, data on the performance characteristics of these tests are limited [9].

Polymerase Chain Reaction (PCR) based methods are sensitive but require expensive thermocycling equipment, electricity, and cold chains for reagents [9–13]. The high sensitivity of nucleic acid amplification technologies (NATs) such as 16S rRNA PCR with a detection limit of 10 fg/μl in *N. gonorrhoeae*-positive clinical samples is attributed to their ability to produce a positive test result from a few copies of DNA or RNA. A real-time amplification assay based on the *N. gonorrhoeae opa* gene had a detection limit of 0.4 bacterial DNA copies [9]. This assay detects asymptomatic *N. gonorrhoeae* infections. However, the performance of this assay has yet to be validated for all *Neisseria* species [9]. Despite their many advantages, the implementation of NAATs by laboratories is expensive, and physical design constraints are a major challenge [9]. In addition, NAAT tests, such as the COBAS AMPLICOR test, produce false-positive results with non-pathogenic *Neisseria* species [14].

The sensitivity of PCR for the detection of *N. gonorrhoeae* varies with different specimens. In a previous study, the sensitivity of PCR for gonorrhoea was 88.3% for urine specimens from men and 94.8% for urethral swab specimens from men [15]. A PCR-based assay specific for the *orf1* gene reported a high sensitivity of 25 fg DNA in clinical samples, which is equivalent to 10 gonococcal cells [14].

Large-scale clinical use of the loop-mediated isothermal amplification (LAMP) test, especially in low- and middle-income countries, has not been achieved. This is due to affordability and challenges with field deployment, including the need for cold chain storage of temperature-sensitive reagents [10,16,17].

In this research, we developed *N. gonorrhoeae* primer sets based on a novel method based on a *de novo* genome mining strategy that identifies repeat sequences in *N. gonorrhoeae* genome for use as both PCR and isothermal amplification coupled with lateral flow readout.

## 2. Materials and methods

### 2.1. IMRS genome mining algorithm

As previously described, primers were designed based on the Identical Multi-Repeat Sequence (IMRS) genome mining algorithm [18]. The *Neisseria gonorrhoeae* (*N. gonorrhoeae*) genome (NZ\_CP012028.1) was used as the reference template to identify identical, repetitive sequence substrings of <30 bases that can be used as forward or reverse primers. The NIH Basic Local Alignment Search Tool (BLAST) was used to evaluate the primer pairs to confirm *N. gonorrhoeae* specificity.

### 2.2. *Neisseria gonorrhoeae* genomic DNA

*Neisseria gonorrhoeae* (*N. gonorrhoeae*) DNA was obtained from the American Type Culture Collection (ATCC) (strain FA1090) (ATCC® 700825DQ Batch No. 70035980) at a concentration of  $\geq 1 \times 10^5$  copies/μL. The stock DNA was diluted to a concentration of 100 pg/μL ( $3.5 \times 10^4$  copies/μL) and, after that, serially diluted 100-fold and 10-fold in Tris-EDTA (TE) buffer (Thermo Fisher Scientific, Waltham, Massachusetts, USA) for the *N. gonorrhoeae* IMRS and 16S rRNA PCR, respectively.

### 2.3. 16S rRNA PCR for *N. gonorrhoeae*

As previously described, the 16S rRNA PCR assays were carried out [19]. Briefly, a reaction mixture containing dNTPs, (Thermo Fisher Scientific, Waltham, Massachusetts, USA) (0.2 mM), forward and reverse primers (0.01 mM for each), Taq Hot-Start DNA polymerase (Thermo Fisher Scientific, Waltham, Massachusetts, USA) (1.25 U), genomic template DNA, 1 μL to a final PCR reaction of 25 μL was used. The cycling parameters for *N. gonorrhoeae* 16S rRNA PCR were as follows: 95 °C for 3 min; 30 cycles of: 95 °C for 30 s, 53 °C for 30 s, 72 °C for 30 s; and 72 °C for 5 min and a final hold of 4 °C. Non-template controls were included in all the PCR assays and were also included in the gel electrophoresis.

#### 2.4. *N. gonorrhoeae* IMRS PCR

The *N. gonorrhoeae* identical multi-repeat sequence (*N. gonorrhoeae*-IMRS) assays were carried out as previously described [18] in a reaction mixture containing dNTPs (Thermo Fisher Scientific, Massachusetts, USA) (0.2 mM), forward and reverse primers (0.01 mM for each), Taq Hot-Start DNA polymerase (Thermo Fisher Scientific, Massachusetts, USA) (1.25 U), genomic template DNA, 1 µL to a final PCR reaction of 25 µL. The cycling parameters for *N. gonorrhoeae*-IMRS assay were as follows: 95 °C for 3 min; 30 cycles of: 95 °C for 30 s, 55 °C for 30 s, 72 °C for 30 s; and 72 °C for 5 min and a final hold of 4 °C.

All PCR products were resolved on a 1% agarose gel visualized on a UV gel illuminator system under ethidium bromide staining. Non-template controls were included in all the PCR assays and were also included in the gel electrophoresis.

#### 2.5. Isothermal IMRS amplification assay

The Isothermal (Iso) IMRS amplification was performed in a 25 µL reaction mixture as previously described [20] and consisted of *Bst* 2.0 polymerase (640 U/mL) (New England Biolabs, Massachusetts, USA), with 1 × isothermal amplification buffer, 3.2 µM forward primer, and 1.6 µM reverse primer (Jigsaw Biosolutions, Bengaluru, India) combined with 10 mM dNTPs (Thermo Fisher Scientific, Massachusetts, USA), 0.4 M Betaine (Sigma-Aldrich, Missouri, USA), molecular-grade water and Ficoll (0.4 g/mL) (Sigma-Aldrich, Missouri, USA). Amplification was carried out at 56 °C for 40 min. Amplified products were visualized by gel electrophoresis in a 1% gel. Non-template controls were included in all the PCR assays and were also included in the gel electrophoresis.

#### 2.6. Lower limit of detection

To assess the lower limit of detection (LLOD), for the *N. gonorrhoeae* IMRS PCR assay, genomic DNA was diluted 100 folds from 10<sup>2</sup> pg/µL (3.5 × 10<sup>4</sup> copies/µL) to 10<sup>-6</sup> pg/µL (<1 copies/µL) and 10 folds from 10<sup>2</sup> pg/µL (3.5 × 10<sup>4</sup> copies/µL) to 10<sup>-2</sup> pg/µL (3.5 copies/µL) for the gold standard 16S rRNA PCR. After that, five replicates of each dilution were used for the assays. Amplification products were visualized on a 1% gel. To determine the LLOD of the *N. gonorrhoeae* 16S rRNA and IMRS PCR, probit analysis was performed using the ratio of successful reactions to the total number of reactions performed for each assay.

#### 2.7. Cloning and characterization of *N. gonorrhoeae*-IMRS amplicons

Gene cloning was performed as previously described [20] to confirm the sequences of amplicons obtained from the *N. gonorrhoeae*-IMRS PCR assay. *Neisseria gonorrhoeae* (*N. gonorrhoeae*) gDNA was amplified using Assembly\_IMRS-F (ttccggatgctc-gagtttttcagcaagatGTTGCCCGCCCGGCTCAAAGG) and Assembly\_IMRS-R (agaatattgtagagatcttctagaaagatATAGCGGATTAA-CAAAAATCAGGACAAGGC) primers. The underlined uppercase sequences correspond to the IMRS primers for amplifying the *N. gonorrhoeae* genome, whereas the non-priming overlap lowercase sequence at the 5' - end of the primers sequence corresponds to the homologous sequences in the cloning vector. The resulting ~223 bp amplicons were resolved on a 2% agarose gel to confirm the fragment size and purified using the PureLink™ PCR purification kit (ThermoFisher). The purified amplicon was then ligated into pJET1.2 blunt vector (ThermoFisher) using the NEBuilder® HiFi DNA Assembly kit (NEB) per the manufacturer's instruction. The resulting NEBuilder HiFi DNA Assembly product was transformed into NEB 5αCompetent *E. coli* (#C2987, NEB) following the manufacturer's instructions. Transformed colonies were randomly selected, DNA extracted and Sanger-sequenced using the universal pJET1.2 forward sequencing primer (CGACTCACTATAGGGAGAGCGGC) and pJET1.2 forward sequencing primer (AAGAA-CATCGATTTTCCATGGCAG). The resulting nucleotides were trimmed and analyzed using SnapGene software (GSL Biotech; available at [snapgene.com](https://www.snapgene.com)), aligned to check for similarity or "clonal" differences, and BLAST was used to check for similarity with the *N. gonorrhoeae* genome.

#### 2.8. *N. gonorrhoeae*-lateral flow assay

To generate a visual readout signal of amplicons, an *N. gonorrhoeae*-Lateral Flow Assay (*N. gonorrhoeae*-LFA) was used as previously described [20]. Briefly, 2.5 µL annealing buffer, dNTPs and NaCl 1.75 µL, MgSO<sub>4</sub> 1.2 µL *N. gonorrhoeae* 5' biotinylated forward primer, *N. gonorrhoeae* Reverse primer, *N. gonorrhoeae* 3' FAM labelled probe, *N. gonorrhoeae* Digotexin labelled probe, Ficoll 400 6.25 µL, *N. gonorrhoeae* internal control sequence DNA 2.5 µL, ISO Amp III enzyme mix 2.0 µL, template DNA 5 µL and molecular grade water 1.45 µL in a final master mix volume of 25 µL was used. Incubation of the LFA strip (Milenia Biotec GmbH, Giessen, Germany) was performed at 65 °C for 1 h and 5 µL of reaction mixture was then added to the LFA strip. After that, 2 drops of running buffer was also added.

#### 2.9. *N. gonorrhoeae* IMRS and 16S rRNA PCR real-time PCR assay

The Quant Studio 5 Real-Time PCR System (with Quant Studio Design and Analysis Desktop Software v1.5.) was used as a reference method [21] for the *N. gonorrhoeae*-16S rRNA PCR and the *N. gonorrhoeae*-IMRS-PCR and for determining the sensitivity of the *N. gonorrhoeae* IMRS and 16S rRNA PCR primers for detecting *N. gonorrhoeae*. Serially diluted (starting concentration of 10<sup>3</sup> genome copies/µL serially diluted 100 folds for *N. gonorrhoeae*-IMRS and starting concentration of 10<sup>4</sup> genome copies/µL serially diluted 10-folds for 16S rRNA) genomic DNA were used as the template. The final real-time PCR master mix volume was 10 µL in triplicate and

consisted of the following components: 5  $\mu$ L SYBR Green qPCR Master Mix (Thermo fisher, Massachusetts, USA), 1  $\mu$ L forward and reverse IMRS primer mix, 2.5  $\mu$ L template genomic DNA and 1.5  $\mu$ L molecular grade water. The amplification cycling conditions were 50 °C for 2 min, 95 °C for 10 min, 40 cycles of 95 °C for 15 s and 60 °C for 30 s.

### 2.10. Clinical samples and ethics statement

Vaginal swab samples collected from a cohort of women aged between 19 and 49 years enrolled in an STI screening study at the Kenyatta National Hospital in Nairobi County in 2022 were used in this study. The *N. gonorrhoeae*-16S rRNA conventional PCR assay was used to confirm *N. gonorrhoeae* infections. The use of these samples and the study proposal was approved by the Mount Kenya University Ethical Review Committee (MKU/ERC/1649).

### 2.11. Data analysis

Graphs were plotted with GraphPad Prism version 7.0 (GraphPad Software, San Diego, CA). The mean and standard deviation values were calculated using Excel 2016. To determine the LLOD of *N. gonorrhoeae*-IMRS and *N. gonorrhoeae*-16S rRNA PCR assays (the concentration at which genomic *Neisseria gonorrhoeae* DNA is detected with 95% confidence), probit regression analyses were performed in Excel 2016. Genomic *Chlamydia trachomatis*, *Trichomonas vaginalis* and *Treponema pallidum* DNA were used to determine the specificity of the *N. gonorrhoeae*-IMRS primers for other STIs. Statistical analyses were performed using a *t*-test of GraphPad Prism version 7.0. for two-tailed distribution.  $P < 0.05$  was considered as significant.

## 3. Results

### 3.1. Design and distribution of IMRS primer sets on *Neisseria gonorrhoeae* genome


Sixteen (16) repeats (Table 1) were identified using the IMRS-based genome mining algorithm, which could be used as forward and reverse primers for an amplification assay. These selected primer pairs (NgIMRS-Primer- F 5'- GTTGCCCCGCCCCGCTCAAAGG-3', NgIMRS-Primer-R 5'- ATAGCGGATTAACAAAAATCAGGACAAGGC - 3') were identified on the genome sequence of *Neisseria gonorrhoeae* 35/02 strain (GenBank accession number CP012028, version: CP012028). As depicted in the Circos plot in Fig. 1A, the primers were found to be present at various loci of the sense and antisense strands, allowing them to serve interchangeably as forward or reverse primers.

### 3.2. Amplification of genomic *Neisseria gonorrhoeae* using identical multi-repeat and gold standard 16S rRNA PCR primers


The primers were used to amplify the serially diluted *N. gonorrhoeae* DNA genome. The IMRS primers could detect *N. gonorrhoeae* genomic DNA from a concentration of as low as 1 fg/ $\mu$ L (Fig. 1B), whereas the gold standard 16S rRNA primers (Fig. 1C) could only detect *N. gonorrhoeae* genomic DNA down to a concentration of 1 pg/ $\mu$ L. This demonstrates that the IMRS-based PCR assay has a higher sensitivity compared to the gold standard 16S rRNA PCR assay. Real-time PCR was also performed using serially diluted genomic *N. gonorrhoeae* DNA as a template and *N. gonorrhoeae*-IMRS primers and conventional *N. gonorrhoeae*-16S rRNA primers. The mean Ct values at each respective dilution were used to plot graphs (Fig. 2C for *N. gonorrhoeae*-16S rRNA primers and for *N. gonorrhoeae*-IMRS primers).

**Table 1**  
*N. gonorrhoeae*-IMRS primer targets on the bacterial genome.

	<i>N. gonorrhoeae</i> -IMRS Primer amplification regions
1	277972-278095
2	442143-442266
3	534149-534272
4	659675-659798
5	741065-741188
6	1170515-1170638
7	1350817-1350941
8	1818578-1818701
9	2064281-2064404
10	505504 - 505381 (-)
11	825946 - 825823 (-)
12	924241 - 924118 (-)
13	960240 - 960117 (-)
14	978398 - 978275 (-)
15	1858932 - 1858809 (-)
16	1971573 - 1971450 (-)

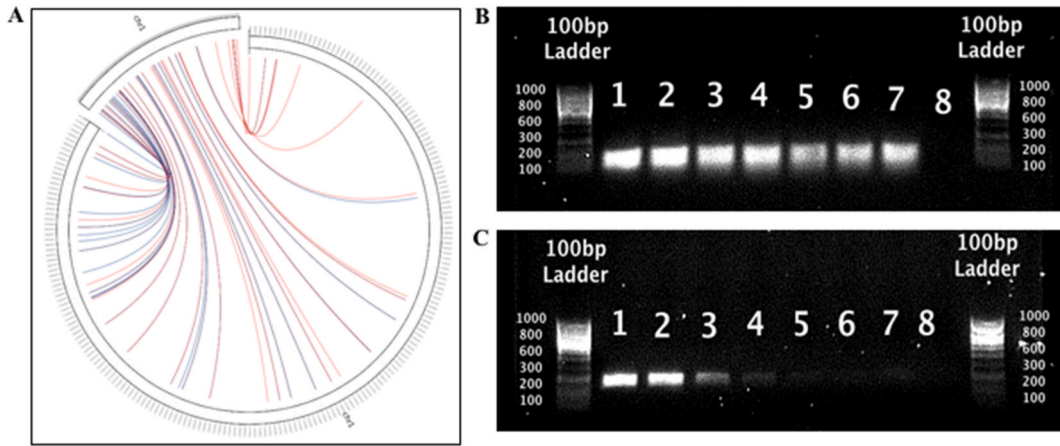

**Fig. 1.** *N. gonorrhoeae*-IMRS Primer targets on the genome and gel image of the *N. gonorrhoeae*-IMRS and 16S rRNA PCR assay. (A) Circos plot for the distribution of identical multi repeat sequence (IMRS) primers in the *Neisseria gonorrhoeae* genome. *Neisseria gonorrhoeae* IMRS primer A (blue lines) and *Neisseria gonorrhoeae* IMRS primer B (red lines) both have 16 repeats. Image of 10-fold serially diluted (1–100, 2–10, 3–1, 4–0.1, 5–0.01, 6–0.01, 7–0.001 and 8 - NTC) (pg/μL) genomic *Neisseria gonorrhoeae* DNA amplicons resolved on a 1% gel using IMRS primers (B) and gold standard 16S rRNA PCR (C). The final length of PCR products obtained using IMRS primers was 124 bp and gold standard 16S rRNA PCR was 260 bp (Supplementary Fig. 1).

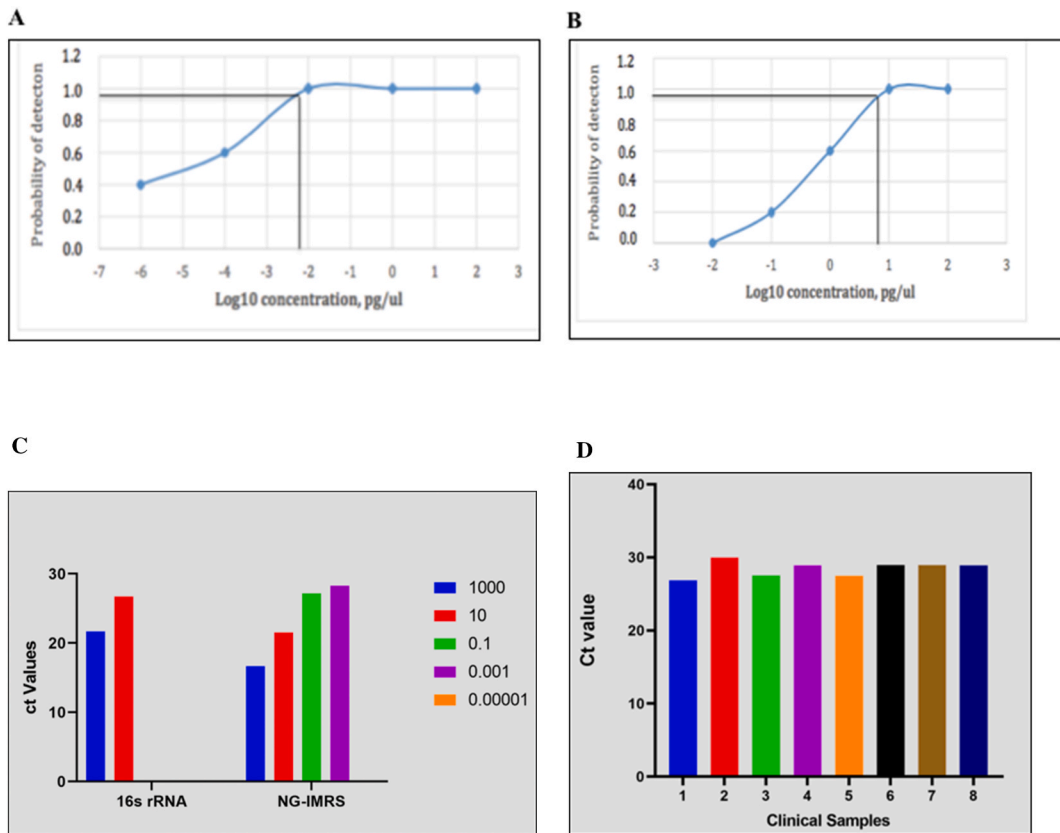

**Fig. 2.** Probit analysis estimation for IMRS. As indicated (A), the IMRS primers for *Neisseria gonorrhoeae* had an LLOD = 6 fg/μL. B: Probit analysis estimation for 16S rRNA PCR. As indicated, gold standard primers for *Neisseria gonorrhoeae* had an LLOD = 4.3096 pg/μL. C: Plots showing Ct values obtained from amplification of genomic *Neisseria gonorrhoeae* serially diluted genomic DNA using RT-PCR for the 16S rRNA and *N. gonorrhoeae*-IMRS primers. The legend shows the starting concentration of the genomic DNA dilutions. D: A total of eight PCR confirmed *Neisseria gonorrhoeae* clinical samples were used to validate the *N. gonorrhoeae*-IMRS primers using RT-PCR.

### 3.3. Calculation of the lower limit of detection (LLOD) for IMRS PCR and gold standard 16S rRNA PCR

To determine the lowest limit of detection (LLOD) of the IMRS PCR relative to the gold standard 16S rRNA PCR, probit statistic was performed using *N. gonorrhoeae* genomic DNA serially diluted 100-fold (Table 2A) and 10-fold (Table 2B) and used as the template for the *N. gonorrhoeae*-IMRS and 16S rRNA PCR respectively. Fig. 2A shows the probit plot for the *N. gonorrhoeae*-IMRS PCR assay, and Fig. 2B shows the probit plot for the gold standard 16S rRNA PCR assay. The LLOD was calculated as the concentration at which *N. gonorrhoeae* DNA can be detected with 95% confidence. Probit analysis estimation for *N. gonorrhoeae*-IMRS PCR, Coefficient  $\chi = -5.0456$ , *P*-value 0.0095 (Table 2C). As indicated, the IMRS primers for *Neisseria gonorrhoeae* had an LLOD = 6 fg/ $\mu$ L, Fig. 2A. Probit analysis estimation for 16S rRNA PCR  $\chi = -7.1656$ , *P*-value 0.9976, Table 2C. As indicated, gold standard primers for *Neisseria gonorrhoeae* had an LLOD = 4.31 pg/ $\mu$ L, Fig. 2B. Table 2C shows the statistics obtained from the Probit analysis indicating that the *N. gonorrhoeae*-IMRS PCR assay had increased sensitivity compared to the gold standard 16S rRNA PCR assay.

### 3.4. Isothermal amplification of cultured *Neisseria gonorrhoeae* DNA

Extracted DNA (38 ng/ $\mu$ L) from cultured *Neisseria gonorrhoeae* was used to perform ISO-IMRS. As shown in Fig. 3A, the reaction products were visualized on a 1% gel.

### 3.5. Validation of *N. gonorrhoeae* IMRS PCR assays with clinical samples

The *N. gonorrhoeae* IMRS primers were validated with DNA isolated from 12 culture confirmed clinical samples. The 16S rRNA PCR assay was run as a control, and amplicons were resolved on a 1% gel, Fig. 3B. As shown in Fig. 3C, the *N. gonorrhoeae* IMRS PCR assay successfully amplified DNA extracted from the 12 culture confirmed samples tested. Using real-time PCR as a reference, three cultured samples were scored as false positive results (Fig. 3D). A total of eight *N. gonorrhoeae* positive DNA samples were used to evaluate the reliability of the *N. gonorrhoeae*-IMRS primers for detecting *N. gonorrhoeae* DNA using RT-PCR assay, as shown in Fig. 2D. The *N. gonorrhoeae*-IMRS primers showed excellent concordance with the results obtained from PCR diagnosis, successfully identifying *N. gonorrhoeae* infections in infected samples.

### 3.6. Cloning and sequencing of *N. gonorrhoeae*-IMRS amplicons


DNA from transformed *E. coli* cells was extracted and sequenced (<https://zenodo.org/records/10224011>). Multiple sequence alignment confirmed that the sequences matched *N. gonorrhoeae* deposited on publicly available resources at *in silico* predicted genomic coordinates. These results confirmed that the *N. gonorrhoeae* IMRS primers were specific for targets within the *N. gonorrhoeae* genome.

As shown in Fig. 4A and B, the *N. gonorrhoeae* -IMRS primers were specific only for *N. gonorrhoeae* DNA. In addition, we designed a probe specific to amplified sequences within the *N. gonorrhoeae* genome. As indicated in Fig. 4C, the *N. gonorrhoeae*-IMRS probe was specific to serially diluted *N. gonorrhoeae* genomic DNA.


**Table 2**

*Neisseria gonorrhoeae* genomic DNA was serially diluted 100 folds (A) and 10-folds (B) and used as template for the *N. gonorrhoeae*-IMRS and 16S rRNA PCR to estimate the lower limit of detection. Table (C) shows the statistics obtained from the Probit analysis.


IMRS PCR		
Serial dilutions (pg/ $\mu$ L)	Replicates (5)	
100		4/5
1		4/5
0.01		4/5
0.0001		3/5
0.000001		2/5
16S rRNA PCR		
Serial dilutions (pg/ $\mu$ L)	Replicates (5)	
100		5/5
10		5/5
1		3/5
0.1		1/5
0.01		1/5
Assay	Coefficient	<i>P</i> -Value
16S rRNA	-7.1656	0.9976
IMRS	-5.0456	0.0095

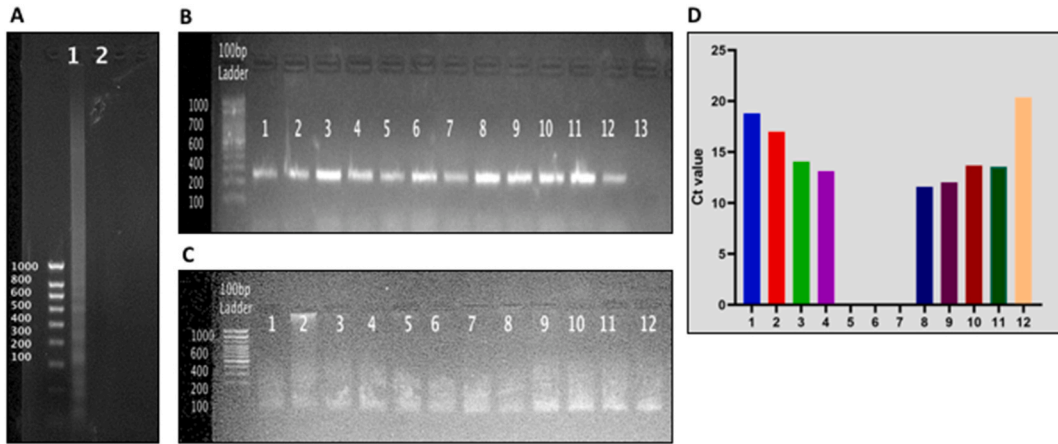

**Fig. 3.** Gel images of DNA amplified using the Isothermal *N. gonorrhoeae*-IMRS primers (A) and the 16S rRNA (B) and the *N. gonorrhoeae* IMRS PCR (C) and Ct values (D) obtained from the RT-PCR assay using cultured *Neisseria gonorrhoeae* clinical isolates. *N. gonorrhoeae*-IMRS primers were used in the RT-PCR. Amplification products were resolved on a 1% gel. 1 – sample, 2, Negative control, (Figure A). 1–12 samples, 13 – Negative control (Figure B and C). The *N. gonorrhoeae*-IMRS assay had a sensitivity of between 66 and 100% and a specificity of 71% (Supplementary Fig. 2).

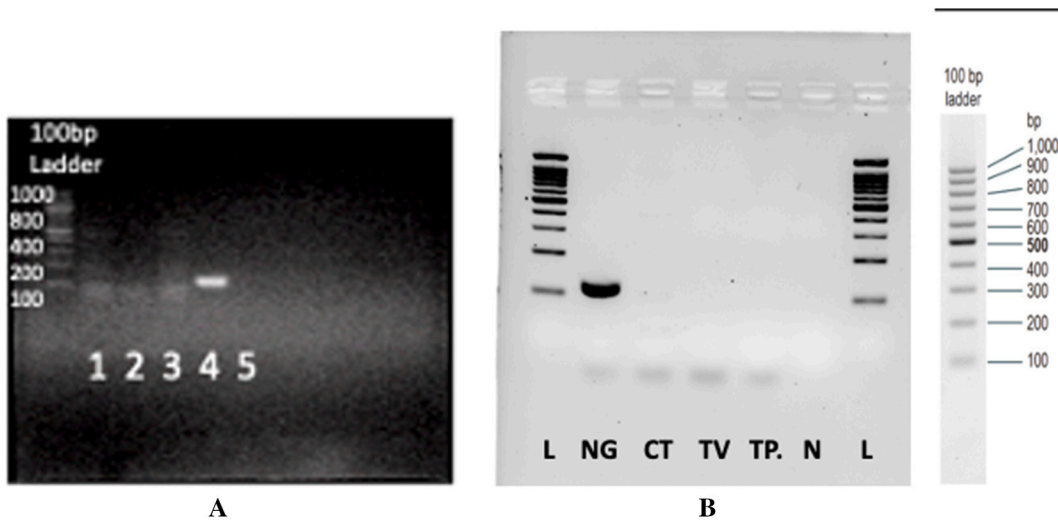

**Fig. 4.** To determine the specificity of *N. gonorrhoeae*-IMRS primers, *Chlamydia trachomatis*, *Treponema pallidum* and *Trichomonas vaginalis* genomic DNA was used as PCR template (A). (B) Gel image of human DNA extracted from venous blood collected from donors (1–3). The extracted DNA served as template for the *N. gonorrhoeae*-IMRS PCR assay and PCR products were resolved on a 1% gel as shown, 4, Positive control – genomic *N. gonorrhoeae* DNA, 5 – Negative control. (C) Ct values from RT-PCR using a probe specific to *Neisseria gonorrhoeae* sequences (Supplementary Fig. 3).

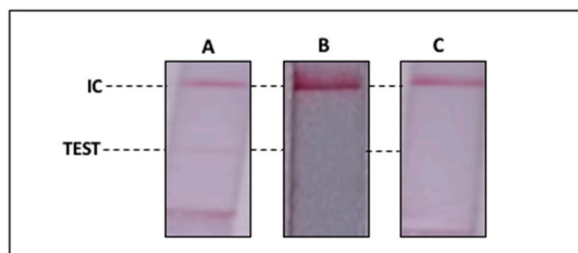

**Fig. 5.** Detection of serially diluted *Neisseria gonorrhoeae* DNA. (A): 100 pg/μL, (B):10 pg/μL, (C): Non Template Control. Strips were incubated at 65 °C for 1 h IC – Internal Control.

### 3.7. *Neisseria gonorrhoeae* lateral flow assay (LFA)

A visual readout of diluted genomic DNA was developed using modified primers and labelled probes, as shown in Fig. 5. The LFA readout of the amplification products successfully indicated the potential for the Iso-IMRS-based assay to be used in field settings.

## 4. Discussion

A key priority in the control of curable STIs is the development of simple and rapid point-of-care tests that help overcome the existing barriers to access and timeliness of diagnosis [22]. In the present study, we used genome-mining approaches that identify identical multi-repeat sequences (IMRS) distributed throughout the *N. gonorrhoeae* genome to design primers that target numerous regions. Compared to the gold standard 16S rRNA PCR, amplifying specific sequences on the *N. gonorrhoeae* genome using IMRS primers was more sensitive, generating many amplicons of varying sizes from the genomic DNA, as shown in Fig. 1B. The IMRS PCR assay was more sensitive with a lower detection limit of 6 fg/ $\mu$ L (2.5 copies/ $\mu$ L, Fig. 2A), representing approximately 1052 times better sensitivity. A previous study based on the Loop-Mediated Isothermal Amplification (LAMP) of Meningococcal genomic DNA reported a lower limit of 10 genome copies per reaction [23]. Compared to the LAMP assay, the *N. gonorrhoeae*-IMRS assay yielded a higher detection limit.

Using real-time PCR to determine the sensitivity of the *N. gonorrhoeae*-IMRS PCR method, we were able to detect up to 1 genome copies/mL of *N. gonorrhoeae* genomic DNA. Efficient asymptomatic *N. gonorrhoeae* infection detection has been identified as a priority in high burden countries. In this perspective, it would be necessary to evaluate the sensitivity of the *N. gonorrhoeae*-IMRS assay in large number of clinical samples.

A study that used real-time PCR assay to target the gonococcal 23S rRNA C2611T and A2059G mutations from vaginal swabs reported Ct value ranges of 29.87–37.98 and 29.2–39.44 respectively [24]. Using clinical samples to validate the *N. gonorrhoeae*-IMRS assay, we obtained Ct value ranges of between 27 and 29 (Fig. 2D), indicating increased sensitivity of the novel assay.

We further showed that isothermal amplification using IMRS (ISO-IMRS) PCR primers detected cultured *N. gonorrhoeae* DNA. The concentration of the initial DNA template was 38 ng/ $\mu$ L (Fig. 3). Our findings are similar to a study that showed a novel ISO-IMRS method was specific for detecting gonococcal strains carrying *penA*-60.001 ceftriaxone-resistant allele [25].

Using *N. gonorrhoeae* genomic DNA as a positive control cross-species specificity for the *N. gonorrhoeae*-IMRS primers was determined using *Chlamydia trachomatis*, *Trichomonas vaginalis* and *Treponema pallidum* genomic DNA (Fig. 4A). Cloning of amplicons followed by sequencing with the plasmid-specific primers yielded the reads aligned within the expected *N. gonorrhoeae* genomic sequences, thus confirming that the amplicons were generated as expected.

In the present study, we have also developed an LFA for diagnosing *N. gonorrhoeae* in low-resource settings. Our *N. gonorrhoeae*-LFA can detect amplicons at a concentration of 100 pg/ $\mu$ L (Fig. 5) visible to the naked eye. The *N. gonorrhoeae*-LFA was more sensitive than a gold nanoparticles colorimetric readout test that detected  $10^6$  CFU/mL in whole live cells in culture [26]. The *N. gonorrhoeae*-LFA can, therefore, replace syndromic management or culture-based clinical methods for diagnosing Gonorrhoeae, particularly in low-resource settings.

## 5. Conclusion

The designed assay can potentially be used to screen asymptomatic infections due to its increased sensitivity and specificity. Our ongoing studies with other STI pathogens indicate that the IMRS algorithm-based genome mining is a platform technology that yields efficient diagnostic biomarkers for STIs and step-change in biomarker discovery in pathogens.

### Funding statement

This research was supported by the Royal Society, Future Leaders African Independent Researchers (FLAIR) Scheme (FLR\R1\201314) to JG.

### Ethical statement

The study proposal was approved by the Mount Kenya University Ethical Review Committee (MKU/ERC/1649).

### Data accessibility

All data generated or analyzed during this study are included in this article [and its supplementary information files]. Underlying raw data has been deposited in Zenodo repository <https://zenodo.org/records/10223536> and <https://zenodo.org/records/10224011>

### Additional information

No additional information is available for this paper.

## CRedit authorship contribution statement

**Clement Shiluli:** Resources, Methodology, Investigation, Formal analysis, Data curation, Conceptualization. **Shwetha Kamath:** Supervision, Investigation, Formal analysis, Data curation, Conceptualization. **Bernard N. Kanoi:** Writing – review & editing, Writing – original draft, Supervision, Formal analysis, Data curation, Conceptualization. **Racheal Kimani:** Writing – review & editing, Software, Resources, Formal analysis, Conceptualization. **Michael Maina:** Writing – original draft, Software, Investigation, Data curation, Conceptualization. **Harrison Waweru:** Supervision, Investigation, Formal analysis, Conceptualization. **Moses Kamita:** Methodology, Investigation, Data curation, Conceptualization. **Ibrahim Ndirangu:** Writing – original draft, Formal analysis, Conceptualization. **Hussein M. Abkallo:** Writing – original draft, Supervision, Software, Resources, Methodology, Investigation, Formal analysis, Conceptualization. **Bernard Oduor:** Writing – original draft, Visualization, Supervision, Resources, Methodology, Investigation, Data curation, Conceptualization. **Nicole Pamme:** Writing – review & editing, Writing – original draft, Supervision, Software, Resources, Project administration, Methodology, Investigation, Formal analysis, Conceptualization. **Joshua Dupaty:** Writing – original draft, Visualization, Supervision, Investigation, Formal analysis, Data curation, Conceptualization. **Catherine M. Klapperich:** Visualization, Validation, Resources, Methodology, Funding acquisition, Formal analysis, Data curation, Conceptualization. **Srinivasa Raju Lola-battu:** Writing – review & editing, Writing – original draft, Supervision, Software, Resources, Investigation, Formal analysis, Data curation, Conceptualization. **Jesse Gitaka:** Writing – review & editing, Writing – original draft, Validation, Supervision, Software, Resources, Investigation, Funding acquisition, Conceptualization.

## Declaration of competing interest

The authors declare that they have no known competing financial interests or personal relationships that could have appeared to influence the work reported in this paper.

## Acknowledgment

We thank Dr. Karen Muthembwa, Mr. Meshack Juma and all clinicians working at the STI Clinic in Kenyatta National Hospital. An earlier version of this article can be found on medRxiv [27].

## Appendix A. Supplementary data

Supplementary data to this article can be found online at <https://doi.org/10.1016/j.heliyon.2024.e27344>.

## References

- [1] J. Whelan, V. Abbing-Karahagopian, L. Serino, M. Unemo, Gonorrhoea: a systematic review of prevalence reporting globally, *BMC Infect. Dis.* 21 (2021), <https://doi.org/10.1186/s12879-021-06381-4>.
- [2] J. Rowley, S. Vander Hoorn, E. Korenromp, N. Low, M. Unemo, L.J. Abu-Raddad, et al., Chlamydia, gonorrhoea, trichomoniasis and syphilis: global prevalence and incidence estimates, 2016, *Bull. World Health Organ.* 97 (2019) 548–562, <https://doi.org/10.2471/BLT.18.228486>.
- [3] M.L. Voth, R.P. Akbari, Sexually transmitted proctitides, *Clin. Colon Rectal Surg.* 20 (2007) 58–63, <https://doi.org/10.1055/S-2007-970201>.
- [4] K. Manavi, H. Young, D. Clutterbuck, Sensitivity of microscopy for the rapid diagnosis of gonorrhoea in men and women and the role of gonorrhoea serovars, *Int. J. STD AIDS* 14 (2003) 390–394, <https://doi.org/10.1258/095646203765371277>.
- [5] E.H. Koumans, R.E. Johnson, J.S. Knapp, M.E. St Louis, Laboratory testing for Neisseria gonorrhoeae by recently introduced nonculture tests: a performance review with clinical and public health considerations, *Clin Infect Dis an Off Publ Infect Dis Soc Am* 27 (1998) 1171–1180, <https://doi.org/10.1086/514994>.
- [6] L. Mukenge-Tshibaka, M. Alary, C.M. Lowndes, E. Van Dyck, A. Guédou, N. Geraldo, et al., Syndromic versus laboratory-based diagnosis of cervical infections among female sex workers in Benin: implications of nonattendance for return visits, *Sex. Transm. Dis.* 29 (2002) 324–330, <https://doi.org/10.1097/00007435-200206000-00003>.
- [7] A. Herring, R. Ballard, D. Mabey, R.W. Peeling, Evaluation of rapid diagnostic tests: chlamydia and gonorrhoea, *Nat. Rev. Microbiol.* 4 (2006), <https://doi.org/10.1038/NRMICRO1562>.
- [8] OPS/OMS, Laboratory diagnosis of sexually transmitted infections, including human immunodeficiency virus, *Clin. Obstet. Gynecol.* 25 (1982) 555–561.
- [9] L.H. Bachmann, R.E. Johnson, H. Cheng, L.E. Markowitz, J.R. Papp, I.W.H. Edward, Nucleic acid amplification tests for diagnosis of Neisseria gonorrhoeae oropharyngeal infections, *J. Clin. Microbiol.* 47 (2009) 902, <https://doi.org/10.1128/JCM.01581-08>.
- [10] L.K. Ng, I.E. Martin, The laboratory diagnosis of Neisseria gonorrhoeae, *Can. J. Infect Dis. Med. Microbiol.* 16 (2005) 15–25, <https://doi.org/10.1155/2005/323082>.
- [11] J.B. Mahony, K.E. Luinstra, M. Tyndall, J.W. Sellors, J. Krepel, M. Chernesky, Multiplex PCR for detection of Chlamydia trachomatis and Neisseria gonorrhoeae in genitourinary specimens, *J. Clin. Microbiol.* 33 (1995) 3049–3053, <https://doi.org/10.1128/jcm.33.11.3049-3053.1995>.
- [12] D.H. Martin, C. Cammarata, B. Van Der Pol, R.B. Jones, T.C. Quinn, C.A. Gaydos, et al., Multicenter evaluation of AMPLICOR and automated COBAS AMPLICOR CT/NG tests for Neisseria gonorrhoeae, *J. Clin. Microbiol.* 38 (2000) 3544–3549, <https://doi.org/10.1128/JCM.38.10.3544-3549.2000>.
- [13] B. Van Der Pol, D.H. Martin, J. Schachter, T.C. Quinn, C.A. Gaydos, R.B. Jones, et al., Enhancing the specificity of the COBAS AMPLICOR CT/NG test for Neisseria gonorrhoeae by retesting specimens with equivocal results, *J. Clin. Microbiol.* 39 (2001) 3092–3098, <https://doi.org/10.1128/JCM.39.9.3092-3098.2001>.
- [14] C.W.M. Geraats-Peters, M. Brouwers, P.M. Schneeberger, A.G.M. Van Der Zanden, S.M. Bruisten, G. Weers-Pothoff, et al., Specific and sensitive detection of Neisseria gonorrhoeae in clinical specimens by real-time PCR, *J. Clin. Microbiol.* 43 (2005) 5653–5659, <https://doi.org/10.1128/JCM.43.11.5653-5659.2005>.
- [15] K.A. Crotchfelt, L.E. Welsh, D. DeBonville, M. Rosenstraus, T.C. Quinn, Detection of Neisseria gonorrhoeae and Chlamydia trachomatis in genitourinary specimens from men and women by a coamplification PCR assay, *J. Clin. Microbiol.* 35 (1997) 1536–1540, <https://doi.org/10.1128/jcm.35.6.1536-1540.1997>.
- [16] OHS, *Infection Surveillance* 70 (2018).

- [17] B. Mlynarczyk-Bonikowska, M. Malejczyk, S. Majewski, M. Unemo, Antibiotic resistance and NG-MAST sequence types of *Neisseria gonorrhoeae* isolates in Poland compared to the world, *Postep Dermatologii i Alergol* 35 (2018) 346–551, <https://doi.org/10.5114/ada.2018.79780>.
- [18] L.S. Raju, S. Kamath, M.C. Shetty, S. Satpathi, A.K. Mohanty, S.K. Ghosh, et al., Genome mining-based identification of identical multirepeat sequences in *Plasmodium falciparum* genome for highly sensitive real-time quantitative PCR assay and its application in malaria diagnosis, *J Mol Diagnostics* 21 (2019) 824–838, <https://doi.org/10.1016/j.jmoldx.2019.04.004>.
- [19] R. Verma, S. Sood, M. Bala, N. Mahajan, A. Kapil, V.K. Sharma, et al., Evaluation of an opa gene-based nucleic acid amplification test for detection of *Neisseria gonorrhoeae* in urogenital samples in North India, *Epidemiol. Infect.* 140 (2012) 2110–2116, <https://doi.org/10.1017/S0950268811002883>.
- [20] N. Kolluri, S. Kamath, P. Lally, M. Zanna, J. Galagan, J. Gitaka, et al., Development and clinical validation of iso-IMRS: a novel diagnostic assay for *P. Falciparum* malaria, *Anal. Chem.* 93 (2021) 2097–2105.
- [21] M. Laig, C. Fekete, N. Majumdar, Digital PCR and the QuantStudio™ 3D digital PCR system, *Methods Mol. Biol.* 2065 (2020) 209–231, [https://doi.org/10.1007/978-1-4939-9833-3\\_16](https://doi.org/10.1007/978-1-4939-9833-3_16).
- [22] J.E. Aledort, A. Ronald, M.E. Rafael, F. Girosi, P. Vickerman, S.M. Le Blancq, et al., Reducing the burden of sexually transmitted infections in resource-limited settings: the role of improved diagnostics, *Nature* 444 (Suppl 1) (2006) 59–72, <https://doi.org/10.1038/NATURE05447>.
- [23] D.K. Lee, E.J. Kim, P.E. Kilgore, S.A. Kim, H. Takahashi, M. Ohnishi, et al., Clinical evaluation of a loop-mediated isothermal amplification (LAMP) assay for rapid detection of *Neisseria meningitidis* in cerebrospinal fluid, *PLoS One* 10 (2015) e0122922, <https://doi.org/10.1371/JOURNAL.PONE.0122922>.
- [24] E. Trembizki, C. Buckley, B. Donovan, M. Chen, R. Guy, J. Kaldor, et al., Direct real-time PCR-based detection of *Neisseria gonorrhoeae* 23S rRNA mutations associated with azithromycin resistance, *J. Antimicrob. Chemother.* 70 (2015) 3244–3249, <https://doi.org/10.1093/JAC/DKV274>.
- [25] K. Shimuta, S. Ichi Nakayama, H. Takahashi, M. Ohnishi, A loop-mediated isothermal amplification assay targeting *Neisseria gonorrhoeae* penA-60.001, *Antimicrob. Agents Chemother.* 64 (2020) 1–19, <https://doi.org/10.1128/AAC.01663-19>.
- [26] T. Oeschger, L. Kret, D. Erickson, Lateral flow assay for detection and recovery of live cell *Neisseria gonorrhoeae*. <https://doi.org/10.1016/j.crbiot.2022.08.002>, 2022.
- [[27]] C. Likhovole, S. Srinivas, B. Kanoi, R. Kimani, M. Maina, H. Waweru, J. Gitaka, in: Improving Gonorrhoea Molecular Diagnostics: Genome Mining-Based Identification of Identical Multi-Repeat Sequences (IMRS) in *Neisseria Gonorrhoeae* Genome, medRxiv, 2023, <https://doi.org/10.1101/2023.07.20.23292966>, 2023-07.